

Absorbent Table Covers

Allesorb® Impervious

CODE	DESCRIPTION (LxW)	CASE QTY	CASE DIMENSIONS
99-UPIMP4030	Allesorb Impervious Table Cover, 40x30" (102x76cm)	10/Bag, 60/CS	18x12x11" (46x30x27cm)
99-UPIMP4060	Allesorb Impervious Table Cover, 40x60" (102x152cm)	1/Pack, 30/CS	15x14x11" (37x36x29cm)
99-UPIMP4090	Allesorb Impervious Table Cover, 40x90" (102x229cm)	1/Pack, 20/CS	16x15x11" (40x37x29cm)
99-UPIMP40120	Allesorb Impervious Table Cover, 40x120" (102x305cm)	1/Pack, 26/CS	18x11x23" (45x29x58cm)
99-UPIMP4030G	Allesorb Premium Impervious Table Cover, 40x30" (102x76cm)	10/Bag, 60/CS	18x12x11" (46x30x27cm)
99-UPIMP4060G	Allesorb Premium Impervious Table Cover, 40x60" (102x152cm)	1/Pack, 30/CS	15x14x11" (37x36x29cm)
99-UPIMP4090G	Allesorb Premium Impervious Table Cover, 40x90" (102x229cm)	1/Pack, 20/CS	16x15x11" (40x37x29cm)
99-UPIMP4090E	Allesorb Impervious Table Cover w/ Elastic, 40x90" (102x229cm)	14	16x15x11" (40x37x29cm)

Materials

- Impervious Polyethylene (PE) Film and Hydrophilic Spunbond Polypropylene (PP) Bi-laminate
- Super Absorbent Polymer Core and Premium Wood Pulp

Product Attributes

- Top layer wicks fluid from surface
- Absorbent middle layer traps and contains fluid
- Impervious base layer prevents strike through
- Allesorb Premium features base layer with improved grip
- Impervious Table Cover w/ Elastic Features a Sewn Elastic Perimeter

Sterilization

Non-Sterile: If a product is delivered non-sterile and is ultimately for use in sterile environments, it is intended to be sterilized by the re-packer. This includes further packaging and sterilization according to the validated processes of the re-packer. These products are to be sterilized by ethylene oxide (EO). Other methods of sterilization have not been validated by Alleset.

Shelf Life

Non-Sterile: Established at 5 years based upon sterile packaged product aging studies and normal warehouse storage conditions. After 5 years, the product should be performance tested to ensure compliance to specifications.

Product Certifications/Standards


- CE certified as class I, according to MDR 2017/745, Annex VIII, Rule 1

Manufacturing Certifications/Registrations

- EN ISO 13485:2016

Storage

- General warehouse conditions

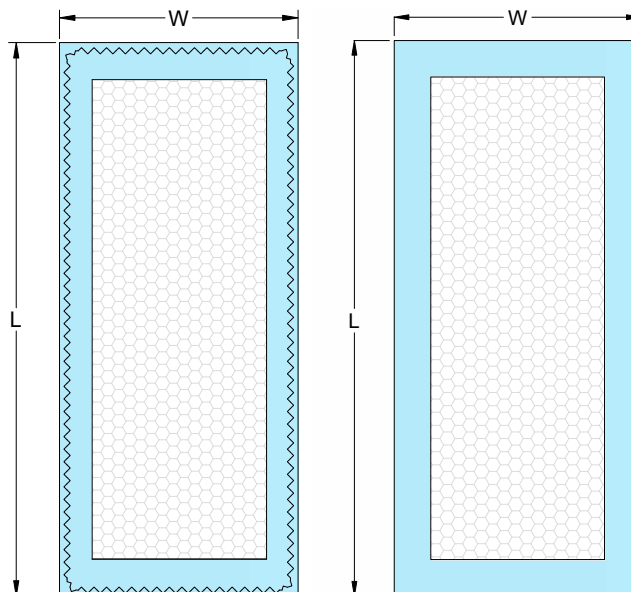

Alleset
4142 Industry Way
Flowery Branch, GA 30542,
USA
www.alleset.com
SM-TDS-002 2022
Country of Origin Abbreviations per ISO 3166-1, Alpha-2 Codes.



GRI Medical & Electronic Technology Co., LTD
1805 Honggao Road,
314031 Jiaxing, Zhejiang
PEOPLE'S REPUBLIC OF CHINA
www.alleset.com
US FDA Establishment Registration #3004911384



GRI-Alleset Limited B.V.
Kingsfordweg 151,
1043 GR Amsterdam,
The Netherlands
www.alleset.com

