

CORE Room Turnaround Kits®

Standard Turnaround Kit

CODE	DESCRIPTION	CASE QTY	CASE DIMENSIONS
49-RTK02	Operating Room (OR) Turnaround Kits	36	24x18x21" (60x45x54cm)

Materials

Absorbent Table Cover

- Impervious Polyethylene (PE) Laminate
- Super Absorbent Polymer Core

Lift Sheet

- Polyethylene Terephthalate (PET)

Armboard Cover

- Polypropylene (PP)/Polyethylene (PE) Laminated Fabric
- Foam Strap with Hook and Loop Closure

Head Rest Cover

- Polypropylene (PP)/Polyethylene (PE) Laminated Fabric

Product Attributes

Absorbent Table Cover:

- Top layer wicks fluid from surface
- Absorbent middle layer traps and contains fluid
- Impervious base layer prevents strike through

Lift Sheet:

- Non-woven PET spunlace construction

Armboard Cover:

- Impervious laminated fabric
- Foam armboard strap included
- Sonic bonded seams

Head Rest Cover:

- Impervious laminated fabric
- Sonic bonded seams

Sterilization

Non-Sterile: If a product is delivered non-sterile and is ultimately for use in sterile environments, it is intended to be sterilized by the re-packer. This includes further packaging and sterilization according to the validated processes of the re-packer. These products are to be sterilized by ethylene oxide (EO). Other methods of sterilization have not been validated by Alleset.

Shelf Life

Non-Sterile: Established at 5 years based upon sterile packaged product aging studies and normal warehouse storage conditions. After 5 years, the product should be performance tested to ensure compliance to specifications.

Product Certifications

- CE certified as Class I, according to MDR 2017/745 Annex VIII, Rule 1

Manufacturing Certifications/Registrations


- EN ISO 13485:2016

Storage

- General warehouse conditions

Kit Components


Component	Product Code	Dimensions	Qty
Absorbent Table Cover	99-UPIMP4090	40x90" (102x229cm)	1
Lift Sheet	50-LFT001	34x66" (86x168cm)	1
Armboard Cover	50-AMBC01-1	13x30" (33x76cm)	2
Head Rest Cover	50-HRC001	14.5x15" (37x38cm)	1

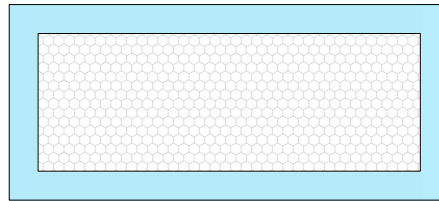


Absorbent Table Cover
99-UPIMP4090



Lift Sheet
50-LFT001

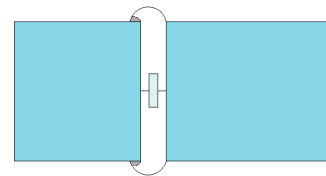

Armboard Cover
50-AMBC01-1

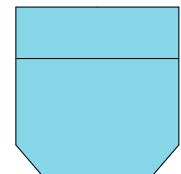

Head Rest Cover
50-HRC001



Alleset
4142 Industry Way
Flowery Branch, GA 30542,
USA
www.alleaset.com
SM-TDS-002 2022
Country of Origin Abbreviations per ISO 3166-1, Alpha-2 Codes.



GRI Medical & Electronic Technology Co., LTD
1805 Honggao Road,
314031 Jiaxing, Zhejiang
PEOPLE'S REPUBLIC OF CHINA
www.alleaset.com
US FDA Establishment Registration #3004911384



GRI-Alleset Limited B.V.
Kingsfordweg 151,
1043 GR Amsterdam,
The Netherlands
www.alleaset.com