

# CORE Room Turnaround Kits®

## Standard Turnaround Kit

CODE	DESCRIPTION	CASE QTY	CASE DIMENSIONS
49-RTK03	Operating Room (OR) Turnaround Kits	36	24x16x17" (60x40x42cm)

### Materials

#### Absorbent Table Cover

- Impervious Polyethylene (PE) Film and Hydrophilic Spunbond Polypropylene (PP) Bi-laminate with Antistatic Treatment
- Super Absorbent Polymer Core

#### Armboard Cover w/ Straps

- Polypropylene (PP)/Polyethylene (PE) Laminated Fabric
- Foam Strap with Hook and Loop Closure

#### Headrest Cover

- Polypropylene (PP)/Polyethylene (PE) Laminated Fabric

### Product Attributes

#### Absorbent Table Cover:

- Top layer wicks fluid from surface
- Absorbent middle layer traps and contains fluid
- Impervious base layer prevents strike through

#### Armboard Cover:

- Impervious laminated fabric
- Foam Armboard Strap Included
- Sonic Bonded seams help prevent strike through

#### Headrest Cover

- Impervious laminated fabric
- Sonic bonded seams

### Sterilization

**Non-Sterile:** If a product is delivered non-sterile and is ultimately for use in sterile environments, it is intended to be sterilized by the re-packager. This includes further packaging and sterilization according to the validated processes of the re-packager. These products are to be sterilized by ethylene oxide (EO). Other methods of sterilization have not been validated by Alleset.

### Shelf Life

**Non-Sterile:** Established at 5 years based upon sterile packaged product aging studies and normal warehouse storage conditions. After 5 years, the product should be performance tested to ensure compliance to specifications.

### Product Certifications/Standards

- CE certified as Class I, according to MDR 2017/745 Annex VIII, Rule 1

### Manufacturing Certifications/Registrations


- EN ISO 13485:2016

### Storage


- General warehouse conditions

### Kit Components


Component	Product Code	Dimensions	Qty
Absorbent Table Cover	99-UPIMP4090	40x90" (102x229cm)	1
Armboard Cover w/ Straps	50-AMBC01	13x30" (33x76cm)	2
Headrest Cover	50-HRC001	14x15" (35.5x38cm)	1

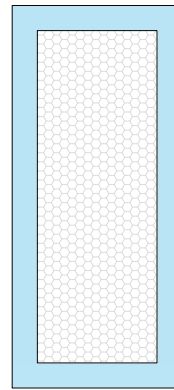

Absorbent Table Cover  
99-UPIMP4090

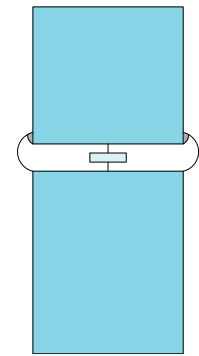

Armboard Cover  
50-AMBC01

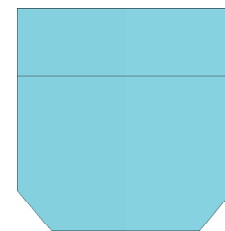

Headrest Cover  
50-HRC001



Alleset  
4142 Industry Way  
Flowery Branch, GA 30542,  
USA  
www.alleset.com  
SM-TDS-003 2023  
Country of Origin Abbreviations per ISO 3166-1, Alpha-2 Codes.



GRI Medical & Electronic Technology Co., LTD  
1805 Honggao Road,  
314031 Jiaxing, Zhejiang  
PEOPLE'S REPUBLIC OF CHINA  
www.alleset.com  
US FDA Establishment Registration #3004911384



Emergo Europe B.V.  
Westervoortsedijk 60,  
6827 AT Arnhem,  
The Netherlands  
www.alleset.com

