

Laparoscopic Accessories

Veress Needle

CODE	DESCRIPTION	CASE QTY	CASE DIMENSIONS
N 37-2150	Veress Needle, 150mm, 14 gauge	500	21x14x13" (53x35x34cm)
S 37-2150-S	Veress Needle, 150mm, 14 gauge	72	6x12x8" (15x31x20cm)
N 37-2120	Veress Needle, 120mm, 14 gauge	500	21x14x13" (53x35x34cm)
S 37-2120-S	Veress Needle, 120mm, 14 gauge	72	6x12x8" (15x31x20cm)

Materials

- Needle: 304 Stainless Steel
- Acrylonitrile Butadiene Styrene (ABS)
- Polycarbonate (PC)
- Polyethylene (PE)

Product Attributes

- **Beveled Edge of Needle:** Helps facilitate needle entry into the abdominal wall.
- **Red-Safety Indicator:** Extends above handle to indicate when the needle tip is exposed.
- **Blunt Tip:** Reduces risk of trauma and helps prevent airflow obstruction.

Sterilization

Sterile: Sterilized by EO (ethylene oxide). Product is sterile unless the package has been damaged or opened.

Non-Sterile: If a product is delivered non-sterile and is ultimately for use in sterile environments, it is intended to be sterilized by the re-packer. This includes further packaging and sterilization according to the validated processes of the re-packer. These products are to be sterilized by ethylene oxide (EO).

Shelf Life

Sterile: 5 years

Non-Sterile: Established at 5 years based upon sterile packaged product aging studies and normal warehouse storage conditions. After 5 years, the product should be performance tested to ensure compliance to specifications.

Product Certifications/Standards


- CE certified as class IIa, according to MDR 2017/745 Annex VIII, Rule 7

Manufacturing Certifications/Registrations

- EN ISO 13485:2016
- Conforms to ISO 594-2:1998: Luer tapers Part 2: Lock fittings

Storage

- General warehouse conditions

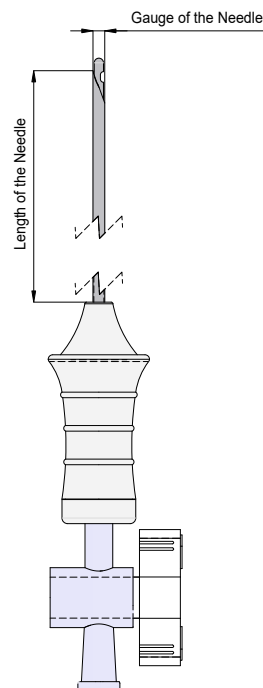

Alleset
4142 Industry Way
Flowery Branch, GA 30542,
USA
www.alleset.com
SM-TDS-003 2023
Country of Origin Abbreviations per ISO 3166-1, Alpha-2 Codes.

GRI Medical & Electronic Technology Co., LTD
1805 Honggao Road,
314031 Jiaxing, Zhejiang
PEOPLE'S REPUBLIC OF CHINA
www.alleset.com
US FDA Establishment Registration #3004911384



Emergo Europe B.V.
Westervoortsedijk 60,
6827 AT Arnhem,
The Netherlands
www.emergobyul.com

